

PRECISA

SEMI-AUTOMATIC TRU-CUT BIOPSY DEVICE WITH TWO RINGS GRIP

The **PRECISA** needle is a device that allows two advancements: 1cm and 2cm.


It allows to perform multiple biopsies by removing the external cannula.

The biopsy under CT doesn't need the introducer: the external cannula (enclosed in the package) in which you insert the stylet with pyramidal tip, is placed into the target.

The safety block avoids accidental cannula advancements. The different steps and control of advancement are performed safely even at different speed thanks to the safety block.

The stereotaxis support relies on the guide allowing to perform multiple biopsies and the use of the safety block.

- Echogenic marker
- Depth marks
- Extra stylet with pyramidal tip
- A two rings grip
- Stereotaxis spacer

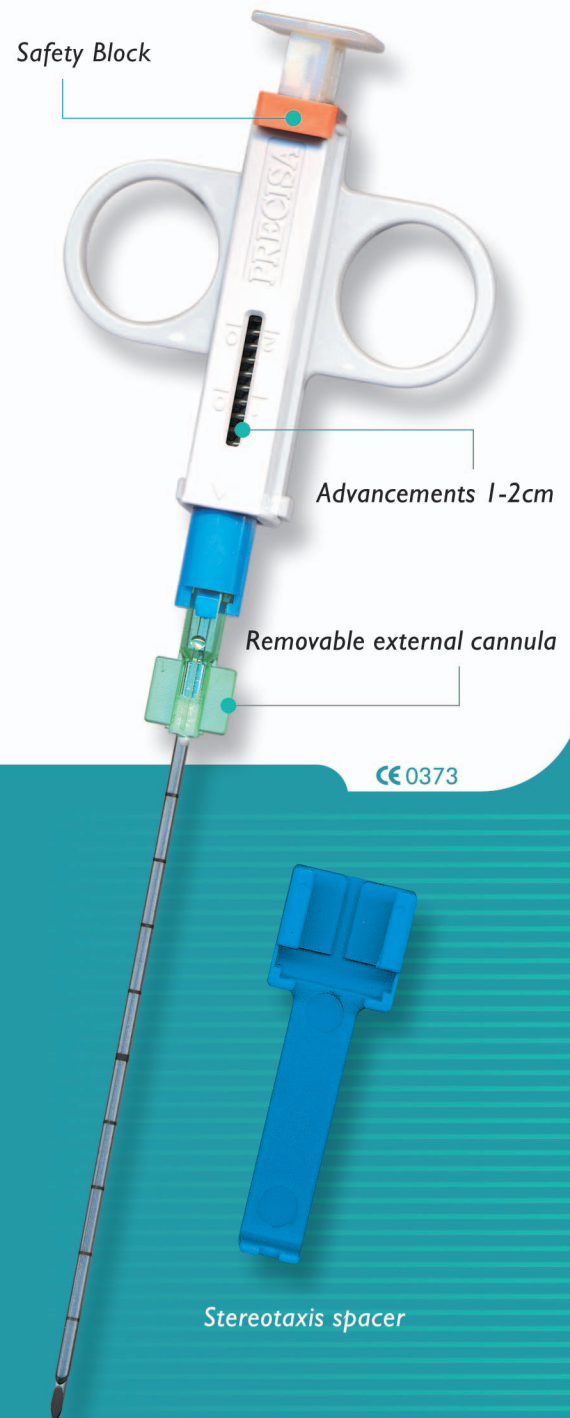

CODE	Diameter	Length	Box
PRE 1310	13G	100mm	5
PRE1410	14G	100mm	5
PRE1415	14G	150mm	5
PRE1420	14G	200mm	5
PRE1610	16G	100mm	5
PRE1615	16G	150mm	5
PRE1620	16G	200mm	5
PRE1625	16G	250mm	5
PRE1807	18G	70mm	5
PRE1810*	18G	100mm	5
PRE1815*	18G	150mm	5
PRE1820	18G	200mm	5
PRE1825	18G	250mm	5
PRE2010*	20G	100mm	5
PRE2015*	20G	150mm	5
PRE2020	20G	200mm	5

* Extra stylet with pyramidal tip

Patent Pending



H.S. HOSPITAL SERVICE S.p.A.
Via delle Valli snc - 04011 Aprilia (LT)
Tel. 069201981 - Fax 0692727871 - Orders fax 069275519
www.hshospitalservice.com e-mail: hs@hshospitalservice.com